

Comparison of phenotypic and genotypic tropism determination in triple-class-experienced HIV patients eligible for maraviroc treatment

Linus Vandekerckhove^{1,2*}, Chris Verhofstede², Els Demecheleer², Stephane De Wit³, Eric Florence⁴, Katrien Fransen⁵, Michel Moutschen⁶, Wendy Mostmans⁷, Kabamba Kabeya³, Nicola Mackie⁸, Jean Plum², Dolores Vaira⁶, Kurt Van Baelen⁷, Ina Vandenbroucke⁷, Veerle Van Eygen⁷, Herwig Van Marck⁷, Dirk Vogelaers¹, Anna Maria Geretti⁹ and Lieven J Stuyver⁷

¹Department of General Internal Medicine, Infectious Diseases and Psychosomatic Medicine, University Hospital Ghent, Ghent, Belgium;

²AIDS Reference Laboratory, University Hospital Ghent, Ghent, Belgium; ³Centre Hospitalier Universitaire Saint-Pierre, Brussels, Belgium;

⁴Department of Clinical Sciences, Institute of Tropical Medicine, Antwerp, Belgium; ⁵AIDS Reference Laboratory, Institute of Tropical

Medicine, Antwerp, Belgium; ⁶Immunology and Infectious Diseases Unit, GIGA-Research, University of Liege, Liege-Sart Tilman, Belgium;

⁷Virco BVBA, Beerse, Belgium; ⁸Imperial College Healthcare National Health Service Trust, London, UK; ⁹University College Medical School and Royal Free Hampstead National Health Service Trust, London, UK

*Corresponding author. Tel: +32-93323398; Fax: +32-933233895; E-mail: linus.vandekerckhove@ugent.be

Received 12 July 2010; returned 20 August 2010; revised 23 October 2010; accepted 30 October 2010

Background: Determination of HIV-1 tropism is a pre-requisite to the use of CCR5 antagonists. This study evaluated the potential of population genotypic tropism tests (GTTs) in clinical practice, and the correlation with phenotypic tropism tests (PTTs) in patients accessing routine HIV care.

Methods: Forty-nine consecutive plasma samples for which an original TrofileTM assay was performed were obtained from triple-class-experienced patients in need of a therapy change. Viral tropism was defined as the consensus of three or more tropism calls obtained from the combination of two independent population PTT assays (Trofile Biosciences, San Francisco, CA, USA, and Virco, Beerse, Belgium), population GTTs and GTTs based on ultra-deep sequencing. If no consensus was reached, a clonal PTT was performed in order to finalize the tropism call. This two-step approach allowed the definition of a reference tropism call.

Results: According to the reference tropism result, 35/49 samples were CCR5 tropic (R5) (patients eligible for maraviroc treatment) and 14/49 were assigned as non-R5 tropic. The non-R5 samples [patients not eligible for maraviroc treatment according to the FDA/European Medicines Agency (EMA) label] group included both the CXCR4 (X4) samples and the dual and mixed CCR5/CXCR4 (R5/X4) samples. Compared with TrofileTM population PTTs, population GTTs showed a higher sensitivity (97%) and a higher negative predictive value (91%), but almost equal specificity and an equal positive predictive value.

Conclusions: In line with recent reports from clinical trial data, our data support the use of population genotypic tropism testing as a tool for tropism determination before the start of maraviroc.

Keywords: HIV-1 tropism, genotypic tests, phenotypic tests

Introduction

Although more than 25 antiretroviral drugs within six different classes have been approved by the US Food and Drug Administration (FDA) and the European Medicines Agency (EMA), current treatment options against HIV-1 remain compromised by broad drug cross-resistance and short- and long-term toxicity.^{1–3} The finding that several commonly prescribed antiretrovirals, including didanosine, abacavir and ritonavir-boosted protease inhibitors, may be associated with increased

cardiovascular risk underscores the need for drugs with an improved safety profile.^{4,5} Entry of HIV-1 into lymphocytes and monocytes requires binding of the envelope glycoprotein (gp120) to the CD4 receptor, followed by interaction with one of two main co-receptors, CCR5 or CXCR4. The preferential use of CCR5 or CXCR4, called viral tropism, is mainly determined by the amino acid sequence of the V3 region of gp120, although regions outside V3, such as V1/V2, C4 and the bridging sheet, may also be involved.⁶ Strains that use the CCR5 co-receptor are called CCR5-tropic or R5 viruses, those that use CXCR4 are

labelled as CXCR4-tropic or X4 viruses. Dually tropic strains can use both CCR5 and CXCR4 and are known as R5/X4 viruses. As the CD4 cell count starts to drop in patients with late stage HIV-1 infection, the prevalence of X4 or R5/X4 strains rises within the viral quasi-species, and starting from a minority population, these strains can finally emerge as a majority population.⁷

Maraviroc is the first CCR5 antagonist to receive FDA and EMEA approval. Determination of viral tropism is a prerequisite to the use of these CCR5 antagonists. Several techniques to determine HIV-1 tropism have been developed over the years, including phenotypic and genotypic methods, which show varying degrees of correlation in direct comparisons. A diagnostic gold standard is currently lacking. Among the phenotypic tropism tests (PTTs), the Trofile™ assay (Monogram Biosciences, San Francisco, CA, USA) in its original version was used for the prospective recruitment of patients within clinical trials.⁸ The original Trofile™ assay (OTA) allowed detection of CXCR4-using strains with a sensitivity approaching 100% at a prevalence of $\geq 10\%$ of the quasi-species in clonal mixtures.⁹ A more sensitive assay version, the 'enhanced sensitivity Trofile™ assay' (ESTA) increased the sensitivity of detection down to 0.3%, again based on mixtures of viral clones.¹⁰ Another PTT has been developed by Virco (Beerse, Belgium) for research purposes only.¹¹ This assay was able to detect CXCR4-using strains in clinical samples when present at $\geq 10\%$ of the quasi-species.¹² It should be noted, however, that the optimal sensitivity of detection of CXCR4-using virus for effectively guiding clinical practice has not been conclusively determined.

Co-receptor use can also be predicted based on the amino acid sequence of the gp120 V3 loop using several publicly available interpretation algorithms.^{13,14}

Among these, Geno2Pheno_[coreceptor] uses a support vector machine statistical learning model. Geno2Pheno has been designed by taking mainly genetic and phenotypic information derived from HIV-1 subtype B into consideration. The recent introduction of the 454 technology for ultra-deep sequencing offers a new opportunity for characterizing the viral quasi-species by generating several thousands of sequences from each sample,^{15–17} thus allowing the fine characterization of drug-resistance patterns.^{15,18}

According to the FDA and EMEA labels, maraviroc can be used only in patients with R5 virus. In the registration trials, the OTA was used for PTTs. The OTA, however, is no longer available, as it was replaced by the ESTA in July 2008. Current consensus from the European, British and German–Austrian guidelines for tropism testing indicate that tropism determination can be performed by a PTT or genotypic tropism test (GTT).^{19–21} The supporting evidence for the GTT as an acceptable alternative to the PTT, and for ESTA as a valid PTT assay, has been mainly derived from retrospective evaluations of patients enrolled in clinical trials.^{10,22,23} The goal of this study was to evaluate the usefulness of an academic developed population GTT in comparison with the OTA (assay available at the time of sampling) in an independent dataset in routine clinical practice.

Methods

Study subjects and clinical samples

Plasma samples were collected from 49 patients eligible for maraviroc initiation, for which an OTA was performed. Patients were selected

from the maraviroc expanded access programme in three Belgian AIDS Reference Centres (Institute for Tropical Medicine, Antwerp; Centre Hospitalier Universitaire Saint-Pierre, Brussels; and University Hospital, Liege), the Royal Free Hospital HIV Centre and Wharfedale Clinic, St Mary's Hospital (London, UK) between 1 October 2007 and 31 December 2008. Patients were triple-class experienced, with advanced HIV disease (median CD4 nadir count 133 cells/mm³) and with a median HIV-1 RNA load of 80 000 copies/mL at the time of tropism testing. Patients were included within the maraviroc expanded access programme (<http://clinicaltrials.gov/ct2/show/NCT00426660>). Optimized background regimens were guided by drug resistance testing. Informed consent was obtained from all patients enrolled, and the study was approved by the Ethics Committee of the University Hospital, Ghent, Belgium. Samples were stored in aliquots at -80°C until further processing. Plasma HIV-1 RNA was quantified with the Amplicor HIV Monitor v1.5 test (Roche Diagnostics Systems, Basel, Switzerland) with a lower limit of detection of 50 RNA copies/mL.

Definition of the reference tropism result (RTR)

The RTR was defined in two steps. In the first step, the RTR criteria were met if the combination of four different assays showed concordance between at least three of the four assays. The four assays comprised: (i) the Virco population PTT (VpPTT), which uses the NH₂-V4 fragment of gp120 on the Antivirogram® platform; (ii) the OTA, which uses the gp160 envelope gene on a PhenoSense platform; (iii) GTT by ultra-deep pyrosequencing on the 454 platform (UDS-GTT) (454 Life Sciences, Branford, CT, USA), combined with interpretation provided by the Geno2Pheno bio-informatic tool (<http://coreceptor.bioinf.mpi-inf.mpg.de/index.php>); and (iv) GTT by population sequencing (pGTT), combined with the Geno2Pheno bio-informatic tool. In the second step, discordant samples that did not reach the criteria of three or more concordant results were analysed using clonal PTT of up to 24 clones in order to allow a final tropism call.

According to the RTR, results were reported as R5 (patients eligible for maraviroc treatment) and non-R5 tropic (patients not eligible for maraviroc treatment according to the FDA/EMEA label). The non-R5 group included the CXCR4 (X4) and the dual and mixed tropic samples.

V3 haplotype

The term V3 haplotype was used to describe a viral population representing individual viruses with an identical V3 sequence fragment.

VpPTT

The assay was performed as described by Van Baelen *et al.*¹¹ Following RNA extraction, the gp120 NH₂-V4 region was amplified in a one-step RT-PCR. The purified NH₂-V4 amplicons were recombined with an HIV-1 backbone deleted for gp120 NH₂-V4. Recombinant plasmids were transfected into 293T cells. The generated virus stocks were used for infection of U87-CD4-CXCR4 and U87-CD4-CCR5 cells. Infection was evaluated by fluorescence readout. Results were reported as R5, non-R5 or failure.

Virco clonal phenotypic tropism testing

The assay was performed essentially as the VpPTT assay, but the starting material consisted of plasmids containing a clonal amplicon of NH₂-V4. Results were reported as R5, X4 or dual tropic (D), or not viable.

Table 1. Overview of tropism assignment according to the two-step approach and concordance of (i) VpPTT; (ii) OTA; (iii) $_{UDS}$ GTT; and (iv) $_{pGTT}$

Tropism assignment step		Samples showing concordance in at least three assays									Samples lacking concordance in three assays					Reference tropism call		
	group = VpPTT/OTA/ _{UDS} GTT/ _{pGTT}	A = R5/R5/R5/R5	B = non-R5/non-R5/non-R5/non-R5	C = R5/non-R5/R5/R5	D = R5/non-R5/non-R5/non-R5	E = non-R5/R5/non-R5/non-R5	F = non-R5/non-R5/non-R5/fail	G = R5/fail/R5/R5	H = non-R5/fail/non-R5/non-R5	I = fail/R5/R5/R5	J = R5/R5/non-R5/fail	K = R5/R5/non-R5/non-R5	L = non-R5/non-R5/R5/R5	M = R5/fail/non-R5/R5	N = non-R5/fail/non-R5/R5	R5	non-R5	
Step 1 (concordant results for ≥ three assays)																		
R5		23		3				3		4						33		
Non-R5			5		1	2	1		2								11	
Step 2 (samples lacking concordance in three assays)																		
R5												1		1			2	
Non-R5											1		1		1			3
Total																35	14	
	Subtype	AI _(n=1) , C _(n=2) , G _(n=1) , FI _(n=1) , B _(n=16) , non-typeable _(n=2)	A _(n=1) , 13cpX _(n=1) , B _(n=3)	B _(n=3)	B _(n=1)	B _(n=1) , 13cpX _(n=1)	B _(n=1)	G _(n=1) , C _(n=2)	AI _(n=1) , 01_AE _(n=1)	G _(n=1) , B _(n=3)	C _(n=1)	B _(n=1)	18_cpX _(n=1)	19_cpX _(n=1)	AI _(n=1)			

Grey shading indicates samples for which at least one assay failed.

OTA

Clinical samples were collected as part of routine patient management and sent to Monogram Biosciences for tropism testing.⁸ Results were reported as R5, non-R5 or failure.

$_{UDS}$ GTT

The assay was performed as described by Vandenbroucke *et al.*²⁴ In summary, viral RNA was extracted from 256 μ L of plasma and reverse transcribed into cDNA, then the V3 loop was amplified in a nested PCR using bar-coded primers (HXB2 positions: forward primer 6986–7012, reverse primer 7520–7540). To maximize the number of input templates and to minimize variation due to PCR drift, seven parallel RT–PCRs were performed and pooled for each sample.²⁵ Bar-coded amplicons were equimolarly pooled and sequenced on the 454 GS-FLX instrument according to the manufacturer's amplicon sequencing protocol (454 Life Sciences, Roche Diagnostics Systems). Hidden Markov models were applied for alignment of the V3 loop. V3 amino acid sequences were

used for tropism prediction. The statistics on the sequence-read numbers from the combined 49 clinical isolates were as follows: average \pm SD = 11 836 \pm 341 (range 2279–21 220). Prediction of the co-receptor tropism was performed using the clonal Geno2Pheno prediction algorithm (<http://coreceptor.bioinf.mpi-inf.mpg.de/index.php>) with a false positive rate (FPR) of 5%.²⁶ The following calls were generated: with $<5\%$ of the virus population labelled X4 (expressed by the total number of sequencing reads), the sample was considered R5; with $>5\%$ of the virus population as a mixed tropic, the sample was considered non-R5.

$_{pGTT}$

RNA was extracted from 200 to 500 μ L of plasma using the High Pure Viral RNA Kit (Roche Diagnostics Systems). Elution was performed with 50 μ L of elution buffer. Nested-PCR amplification of a fragment spanning the V1 to V4 region of the *env* gene was performed (HXB2 primer positions: outer PCR, forward primer 6540–6560 and reverse primer 7701–7721; and inner PCR, forward primer 6561–6580 and reverse primer 7645–7667).

Sequencing reactions were prepared using the BigDye® Terminator Cycle Sequencing kit v3.1 (Applied Biosystems, Foster City, CA, USA) with degenerate internal primers: sense 5'-AGYRCAGTACAATGYACACATGG-3'; sense 5'-TCAACHCAAYTRCTGTAAATGG-3'; and antisense 5'-ATTCTGGRTCYCKCCTG-3'. Sequencing products were run on an ABI3130xl automated sequencer. Sequence editing and contig assembly were performed using the Smartgene™ HIV software package (Integrated Database Network System, Smartgene, Zug, Switzerland). V3 nucleic acid sequences were used for tropism prediction with the clonal Geno2Pheno prediction algorithm using an FPR of 5%.

Results

Results of four tropism assays

Results of the four assays are summarized in Table 1. Samples were divided into subgroups according to the individual test results (in the following order: V_pPTT/OTA/_{UDS}GTT/_pGTT): A = R5/R5/R5/R5; B = non-R5/non-R5/non-R5/non-R5; C = R5/non-R5/R5/R5; D = R5/non-R5/non-R5/non-R5; E = non-R5/R5/non-R5/non-R5; F = non-R5/non-R5/non-R5/fail; G = R5/fail/R5/R5; H = non-R5/fail/non-R5/non-R5; I = fail/R5/R5/R5 (samples were categorized as showing concordance in at least three assays); J = R5/R5/non-R5/fail; K = R5/R5/non-R5/non-R5; L = non-R5/non-R5/R5/R5; M = R5/fail/non-R5/R5; and N = non-R5/fail/non-R5/R5 (samples were categorized as lacking concordance in at least three assays).

Samples showing concordance in at least three assays

In the groups A–I, at least three assay test results could be obtained. For these samples, a detailed analysis was performed of the discordant call. In group C, three samples were non-R5 by OTA but R5 by the other three assays. The ultra-deep sequencing assay did not detect X4 strains in two of the three samples, and detected a low prevalence (1.6%; but because the cut-off of 5% labelled as R5 only) of X4 strains in the third sample. Thus, the final call for group C favoured R5. In group D, one sample was R5 by V_pPTT but non-R5 by the other three tests. The ultra-deep sequencing assay identified 46% X4 strains. Thus, the final call for group D favoured non-R5. In group E, two samples were R5 by OTA but non-R5 by the other three tests. The ultra-deep sequencing assay identified 24% and 65% X4 strains. Thus, the final call for group E strains favoured non-R5.

Samples lacking concordance in at least three assays

With samples lacking concordance in three or more assays, clonal phenotyping was performed to gain further insights into the discrepancies (see Table 2). In group J, one sample was R5 by V_pPTT and OTA, whereas _{UDS}GTT detected 16% X4 strains, while no result was obtained by _pGTT. By _{UDS}GTT, three different populations (present >5%; haplotypes 1–3) could be distinguished. In one population representing 11% X4 virus of the total population by _{UDS}GTT, clonal phenotyping determined three clones as dual tropic, classifying the sample as non-R5. The subgroup K sample was classified R5 by both V_pPTT and OTA, but X4 by _{UDS}GTT (77% X4 strains) and _pGTT. Of the three populations, two (haplotypes 1 and 2; both representing ±36% of the population) were classified X4 by the G2P prediction tool. However, both populations harboured one or

Table 2. Details of the step 2 tropism identification

Group	V _p PTT	OTA	_{UDS} GTT	_p GTT	Number of haplotypes	Haplotype number	Haplotype V3 amino acid sequence	Populations present >5%	V3 haplotype tropism prediction	Clonal phenotype			Final tropism call
										R5	D	NR	
J	R5	R5	non-R5	fail	32	1	CTRPNNTRKR IHLGPGPGRAFATGGIVGDIRKAYC	11.0	X4	0	3	1	non-R5
						2	-----RS--M^--K-----D--I-----	59.4	R5	8	0	3	
						3	-----RS--M^--K-----D--I-----	15.5	R5	2	0	2	
K	R5	R5	non-R5	non-R5	29	1	CTRPNNTRDKIGIGPGQMFYATQNIIGNIREAHC	35.7	X4	1	0	4	R5
						2	-----E-----H-----	36.0	X4	3	0	3	
						3	-----E-----	21.2	R5	1	0	6	
L	non-R5	non-R5	R5	R5	19	1	CTRPNNTRKGTHIGPGQAFYATDQVIGNIRQAHC	14.9	R5	0	0	1	non-R5
						2	-----N--R--R-----S-----KI--D-----	82.3	R5	0	20	2	
						3	CTRPNNTRKRTHIGPGRAFVTTSALIGNIRQAHC	5.2°	X4	0	0	0	R5
M	R5	fail	non-R5	R5	31	1	-----N-----Q--YA-G-----K-----	7.2	R5	2	0	0	
						2	-----N-----Q--YA-G-----N-----	82.4	R5	18	0	2	
						3	CTRPNNTRKFVRI GPGGPGHAFYTNDIIGDIRKACH	9.3	X4	0	1	1	non-R5
N	non-R5	fail	non-R5	R5	23	1	-----I-----G-H-----Q^--L-----Y-----	9.1	R5	0	1	0	
						2	-----I-----S-H-----Q^--L-----Y-----	75.9	R5	0	0	0	
						3	-----I-----S-H-----Q^--L-----Y-----						

D, dual tropic virus; NR, non-reportable.

°population present at the detection limit of the clonal phenotyping.

more R5 clones determined by clonal phenotyping and no X4 or dual tropic clones. These data suggest that the prediction tool is making an error both on a population as well as on a UDS_{GTT} level. The final call for this sample is R5. The group L sample was non-R5 by both V_{pPTT} and OTA, but R5 by both UDS_{GTT} and p_{GTT} . The two major populations representing haplotype 1 (14.9%) and haplotype 2 (82.3%) were classified R5 by G2P both on a UDS_{GTT} and population genotype level. The clonal phenotyping determined 20 clones as dual tropic from the major population representing 82.3%. Again these data suggest that the prediction tool is making an error both on a population as well as on an UDS_{GTT} level, which might be due to the relatively uncommon subtype (18_cpx). The final call for this sample is non-R5. In group M, one sample was R5 by V_{pPTT} and p_{GTT} , but non-R5 (5.2% X4) by UDS_{GTT} . In clonal phenotyping, 20 clones were classified R5, and no X4 clones were detected. The final call for this sample was R5. Of the three major populations, one population representing 5.2% was not picked up by the clonal phenotyping. Although 20 viable clones were analysed, the clonal phenotyping has a detection limit of around 5% and it cannot be ruled out that the 5.2% virus population was therefore not picked up by the random sampling (see Table 2). In group N, one sample was non-R5 by V_{pPTT} and UDS_{GTT} (10.2% X4), but was R5 by p_{GTT} . By UDS_{GTT} , three different major populations (>5%) could be distinguished. In one population (haplotype 1) representing 9.3% X4 virus of the total population by UDS_{GTT} , clonal phenotyping confirmed this population as dual tropic, classifying the sample as non-R5. The final conclusion of the combination of first-step and second-step analysis was that 35 samples were classified as R5 and 14 samples as non-R5 (Table 1).

Parameters of phenotyping and genotyping methods using RTR

By using the RTR, sensitivity and specificity of the V_{pPTT} , OTA, UDS_{GTT} and p_{GTT} could be calculated, as shown in Table 3. The higher failure rate of the OTA affected the negative predictive value. In general, the sensitivity of the OTA, UDS_{GTT} and p_{GTT} was within the same range.

We evaluated the approximate limit of detection for minority variants using p_{GTT} . Only a few samples harboured a population present between 10% and 25%, which is the estimated detection cut-off for minority variants for population genotyping methods. Therefore, it is hard to come up with a well-defined cut-off from this study. If the minority population drops below 10% the p_{GTT} is not picking up this population, as is shown for sample N (X4 present as 9.3% minority), whereas if the population is above 10% the p_{GTT} is picking up the population (group H, sample with 01_AE subtype, X4 present at 12.5% of the total population).

Subsequently the impact of the HIV subtype on the OTA and p_{GTT} was evaluated. In general, discordances for both the OTA and p_{GTT} are evenly distributed between B and non-B subtypes, but the number of patient samples included was too small to draw conclusions (see Table 1). The concordance with the RTR was 31/36 (86%) for OTA and 34/36 (94%) for p_{GTT} . Of the five discordant OTA samples (C and E), four were subtype B. Of the two discordant GTT samples (K and L), one was subtype B.

Clinical decisions based on OTA

In order to better understand the potential value of the different tropism tests, we performed a theoretical evaluation of the impact of using OTA versus p_{GTT} for guiding clinical decisions (see Figure 1). OTA classified 31/49 samples as R5, allowing 31 patients to start a chemokine receptor antagonist as part of their treatment regimen. According to the definition of reference result described previously, 28 of the 31 samples were correctly identified as R5. Assuming use of p_{GTT} instead of OTA, a total of 36 patients would have been returned as R5, leading to maraviroc use, and of these, 34 could be considered as correctly identified as R5. The net result of adopting p_{GTT} as a diagnostic tool would be to allow an additional five patients to receive maraviroc. Conversely, two patients would have started maraviroc in the presence of X4.

According to the reference tropism test, a total of 35 patients would have been eligible for maraviroc therapy. These 35 eligible individuals were distributed in the OTA results as follows: 28 in the R5 cohort, 4 from the OTA failures and 3 from the non-R5 cohort. Hence, a total of seven patients missed the opportunity to be treated with a co-receptor antagonist, while another three were labelled as candidates for such treatment while non-R5 virus was present.

Similarly, the 35 eligible individuals were distributed in the p_{GTT} results as follows: 34 in the R5 cohort and 1 from the non-R5 cohort. This would result in one patient missing the opportunity to be treated with a co-receptor antagonist, while another two were labelled as candidates for such treatment while X4 virus was present (Table 1).

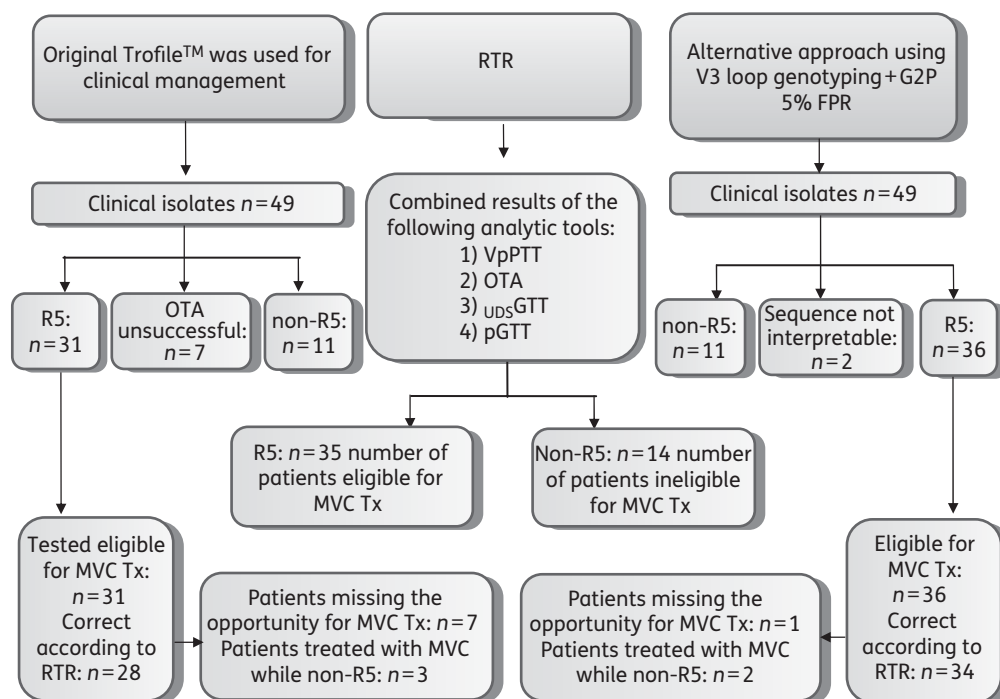
Discussion

In our study, 35 of the 49 plasma samples were classified R5 according to p_{GTT} and 31 samples were classified as R5 according to OTA. If samples for which one assay failed (groups F, G, H, I, J, M and N) were ignored, concordance for all four assays was achieved in 28/36 samples (77%). The concordance with the RTR was 31/36 (86%) for OTA and 34/36 (94%) for p_{GTT} . Discordant results can be attributed to several reasons. In group C, three plasma samples were R5 by p_{GTT} and V_{pPTT} , but non-R5 by OTA. In ultra-deep sequencing, G2P scored 0% of the reads as X4 in two samples. Thus, the FPR of OTA seems to be higher than 0%. Indeed, the only reported FPR was measured using plasma samples of other non-retroviruses.⁸ Another factor that might have driven discordances is that OTA was running as a routine diagnostic test and the other tests were performed retrospectively in a research setting. In other cases, variability of detection at the lower end of assay sensitivity may explain the discrepancy. One sample in group C was R5 by p_{GTT} and V_{pPTT} , but non-R5 by OTA, and ultra-deep sequencing detected 1.6% X4 reads, not accounted for because of the cut-off of 5%. Although it is generally accepted that the detection of minorities is superior for ultra-deep sequencing, the clonal phenotyping was used to generate the final call. By doing this we allowed for the pick up of potential errors in the interpretation tool that would classify both the population and ultra-deep V3 loop sequences identical. In both samples K and L the G2P interpretation algorithm is making a different interpretation than the clonal phenotyping. In these two samples, it is not a question of

Table 3. Characteristics of the genotypic and phenotypic assays

RTR	N	OTA			VpPTT			UDSGTT			pGTT		
		R5	non-R5	F	R5	non-R5	F	R5	non-R5	F	R5	non-R5	F
R5	35	28	3	4	31	0	4	33	2	0	34	1	0
non-5	14	3	8	3	2	12	0	1	13	0	2	10	2
Total	49	31	11	7	33	12	4	34	15	0	36	11	2
PPV		90.3			93.7			97.1			94.4		
NPV		72.7			100			86.7			90.9		
Sensitivity		90.3			100			94.3			97.1		
Specificity		72.7			85.7			92.9			83.3		
Availability		replaced by ESTA			not available			in core facilities			local, home brew protocol		
Regulatory level		CAP/CLIA			research only			research only			dependent on local settings		
TAT		>3 weeks			not relevant			2 weeks			days		
Cost		moderate to high			not relevant			moderate to high			low		

TAT, turnaround time; CAP/CLIA, College of American Pathologists/Clinical Laboratory Improvement Act; PPV, positive predictive value; NPV, negative predictive value; RTR, reference tropism result.

**Figure 1.** Overview of the study. RTR, reference tropism result; MVC, maraviroc; Tx, treatment.

whether or not the ultra-deep sequencing has a better sensitivity for the detection of minorities than the clonal phenotyping, as the populations for which the G2P differs are not present as a minority, but as a majority. The high sensitivity for minority detection of X4 variant by UDSGTT could be of value for the tropism determination of sample M, now classified as R5, as the 5.2% UDSGTT determined X4 could easily be missed by the clonal phenotyping of 24 clones.

A further item arises from the interpretation of sequencing results analysed, especially for rare and complex sequences. In our analysis, if we excluded the samples where one assay failed, the results of the pGTT analysis could not be confirmed by RTR for samples K (subtype B) and L (18_cpx). The results of the OTA could not be confirmed by the RTR for samples C and E (subtype B, $n = 4$; 13_cpx, $n = 1$). Although the ESTA is now performed with new primers that improve tropism prediction

on non-B viral strains, the impact of introducing this new assay with adapted primers might have had a minor influence on the concordance with the RTR, as most samples were subtype B.²⁷

In this study we evaluated the potential value of population genotyping in clinical practice. Importantly, 36 patients could have started and benefited from maraviroc. Detailed analysis suggests that 33 of these 36 patient samples did indeed harbour <5% dual or mixed tropic virus, which is associated with a good response to therapy, shown in a previous analysis by McGovern *et al.*²² Our data support the findings reported earlier from the International HIV Drug Resistance Meeting 2010 indicating that, in treatment-experienced patients, population genotyping is comparable in predicting clinical outcome and response to chemokine receptor antagonist independent of the use of the G2P or position-specific scoring matrix (PSSM) interpretation algorithm with the OTA.²⁸ Our data represent the first dataset in clinical practice where V3 loop tropism determination has been assessed outside the context of a clinical trial with subsequent detailed 454 analysis. In line with the recent European, British and German guidelines, we provide further data to support that population V3 loop genotyping is comparable to OTA in predicting tropism in clinical practice.^{19–21} Population V3 loop genotyping can serve as a tropism determination tool, but full standardization of the sequence method and interpretation remains a requirement. Although our study evaluated two phenotypic assays, population sequencing and ultra-deep V3 loop sequencing, and for some discordant samples also clonal phenotyping, results should be interpreted with caution, as this is a retrospective evaluation of pGTT on a limited number of samples. Prospective evaluation of pGTT in clinical practice remains a topic of further research.

Acknowledgements

We would like to thank Martin Fisher, Margaret Johnson, Nicky Perry and Angela Strang for their contributions.

Funding

This work was supported by internal funding without pharmaceutical sponsoring. L. V. is supported by the Flemish Fund for Scientific Research.

Transparency declarations

A. M. G. receives consultancy fees and research funding from and is on the speaker's bureau for Pfizer, GSK, ViiV, Monogram Biosciences, MSD, Tibotec, J&J, Roche Diagnostics, Gilead, BI, BMS, Virco and Abbott Diagnostics. L. V. receives consultancy fees and research funding from and is on the speaker's bureau for Pfizer, GSK, ViiV, Gilead, BMS, Virco, Tibotec, J&J and MSD. N. M. receives consultancy fees from and is on the speaker's bureau for Tibotec, BMS, ViiV, Merck and Abbott. S. D. W. receives consultancy fees from and is on the speaker's bureau for Pfizer, GSK, BMS and Abbott. E. F. received consultancy fees from Abbott, BMS, Gilead, MSD, J&J and ViiV Healthcare. W. M., K. V. B., I. V., V. V. E., H. V. M. and L. J. S. are employees of Tibotec Virco Virology BVBA. All other authors: none to declare. While Tibotec is commercializing antiviral compounds for HIV-1 treatment, the current work does not represent or influence any of these compound-related commercial activities. Virco has developed and commercialized vicroTYPE for reverse transcriptase and protease resistance testing. Tropism testing tools and

algorithms were developed for research purposes only, and no activities for commercialization are planned or intended.

References

- 1 Llibre JM, Schapiro JM, Clotet B. Clinical implications of genotypic resistance to the newer antiretroviral drugs in HIV-1-infected patients with virological failure. *Clin Infect Dis* 2010; **50**: 872–81.
- 2 Feeney Eoin R, Mallon Patrick WG. Impact of mitochondrial toxicity of HIV-1 antiretroviral drugs on lipodystrophy and metabolic dysregulation. *Curr Pharm Des* 2010; in press.
- 3 Lang S, Mary-Krause M, Cotte L *et al.* Impact of individual antiretroviral drugs on the risk of myocardial infarction in human immunodeficiency virus-infected patients: a case-control study nested within the French Hospital Database on HIV ANRS cohort CO4. *Arch Intern Med* 2010; **170**: 1228–38.
- 4 Friis-Moller N, Reiss P, Sabin CA *et al.* Class of antiretroviral drugs and the risk of myocardial infarction. *N Engl J Med* 2007; **356**: 1723–35.
- 5 Sabin CA, Worm SW, Weber R *et al.* Use of nucleoside reverse transcriptase inhibitors and risk of myocardial infarction in HIV-infected patients enrolled in the D:A:D study: a multi-cohort collaboration. *Lancet* 2008; **371**: 1417–26.
- 6 Cormier EG, Dragic T. The crown and stem of the V3 loop play distinct roles in human immunodeficiency virus type 1 envelope glycoprotein interactions with the CCR5 coreceptor. *J Virol* 2002; **76**: 8953–7.
- 7 Tsiibris AM, Kuritzkes DR. Chemokine antagonists as therapeutics: focus on HIV-1. *Annu Rev Med* 2007; **58**: 445–59.
- 8 Whitcomb JM, Huang W, Fransen S *et al.* Development and characterization of a novel single-cycle recombinant-virus assay to determine human immunodeficiency virus type 1 coreceptor tropism. *Antimicrob Agents Chemother* 2007; **51**: 566–75.
- 9 Petropoulos C, Limoli K, Whitcomb J *et al.* Validation studies defining the performance of Monogram's HIV coreceptor tropism assay. In: *Abstracts of the Fifth European HIV Drug Resistance Workshop, Cascais, Portugal, 2007*. Abstract 57.
- 10 Saag M, Heera J, Goodrich J. Reanalysis of the MERIT study with the enhanced Trofile assay (MERIT-ES). In: *Abstracts of the Forty-eighth Interscience Conference on Antimicrobial Agents and Chemotherapy, Washington, DC, 2008*. Abstract H-1232a. American Society for Microbiology, Washington, DC, USA.
- 11 Van Baelen K, Vandenbroucke I, Rondelez E *et al.* HIV-1 coreceptor usage determination in clinical isolates using clonal and population-based genotypic and phenotypic assays. *J Virol Methods* 2007; **146**: 61–73.
- 12 Vandekerckhove L, Verhofstede C, Vandenbroucke I *et al.* on behalf of the Tropism Study Group. Comparison of population and 454 genotyping and Trofile phenotyping for tropism determination in patients selected for maraviroc initiation. In: *Abstracts of the Eighteenth International HIV Drug Resistance Workshop, Fort Myers, FL, 2009*. Abstract 57.
- 13 Chan SY, Speck RF, Power C *et al.* V3 recombinants indicate a central role for CCR5 as a coreceptor in tissue infection by human immunodeficiency virus type 1. *J Virol* 1999; **73**: 2350–8.
- 14 Jensen MA, van 't Wout AB. Predicting HIV-1 coreceptor usage with sequence analysis. *AIDS Rev* 2003; **5**: 104–12.
- 15 O'Meara D, Wilbe K, Leitner T *et al.* Monitoring resistance to human immunodeficiency virus type 1 protease inhibitors by pyrosequencing. *J Clin Microbiol* 2001; **39**: 464–73.
- 16 Ronaghi M. Pyrosequencing sheds light on DNA sequencing. *Genome Res* 2001; **11**: 3–11.

- 17** Shi MM. Enabling large-scale pharmacogenetic studies by high-throughput mutation detection and genotyping technologies. *Clin Chem* 2001; **47**: 164–72.
- 18** Mitsuya Y, Varghese V, Wang C *et al.* Minority human immunodeficiency virus type 1 variants in antiretroviral-naïve persons with reverse transcriptase codon 215 revertant mutations. *J Virol* 2008; **82**: 10747–55.
- 19** Vandekerckhove LPR, Wensing AMJ, Kaiser R *et al.* Consensus statement of the European guidelines on clinical management of HIV-1 tropism testing. In: *Abstracts of the Tenth International Congress on Drug Therapy in HIV Infection, Glasgow, Scotland, 2010*. Abstract O121.
- 20** Geretti A, Mackie N. *British HIV Association Guidelines on Determining HIV-1 Tropism in Routine Clinical Practice*. <http://www.bhiva.org/Tropism.aspx> (October 2010, date last accessed).
- 21** *German Recommendations for Determining HIV-1 Coreceptor Usage*. <http://www.viro.med.uni-erlangen.de/nrz/recommendation080324.pdf> (September 2010, date last accessed).
- 22** McGovern R, Dong W, Mo T *et al.* Optimization of clinically relevant cutpoints for the determination of HIV co-receptor usage to predict maraviroc responses in treatment experienced (TE) patients using population V3 genotyping. In: *Abstracts of the EACS Conference, Cologne, 2009*. Abstract PE3.4/8.
- 23** Su Z, Gulick RM, Krambrink A *et al.* Response to vicriviroc in treatment-experienced subjects, as determined by an enhanced-sensitivity coreceptor tropism assay: reanalysis of AIDS clinical trials group A5211. *J Infect Dis* 2009; **200**: 1724–8.
- 24** Vandembroucke I, Van Marck H, Mostmans W *et al.* HIV-1 V3 envelope deep sequencing for clinical plasma specimens failing in phenotypic tropism assays. *AIDS Res Ther* 2010; **7**: 4.
- 25** Vandembroucke I, Eygen VV, Rondelez E *et al.* Minor variant detection at different template concentrations in HIV-1 phenotypic and genotypic tropism testing. *Open Virol J* 2008; **2**: 8–14.
- 26** Sing T, Low AJ, Beerenwinkel N *et al.* Predicting HIV coreceptor usage on the basis of genetic and clinical covariates. *Antivir Ther* 2007; **12**: 1097–106.
- 27** Napolitano L, Tressler T, Coakley E *et al.* Incorporation of optimized primers into the Trofile assay substantially improves determination of viral tropism in genetically diverse HIV subtypes. In: *Abstracts of the Forty-ninth Interscience Conference on Antimicrobial Agents and Chemotherapy, San Francisco, 2009*. Abstract 908. American Society for Microbiology, Washington, DC, USA.
- 28** Swenson LC, McGovern AR, Dong W *et al.* Phenotypic screening for HIV tropism versus both population-based and “deep” sequencing. In: *Abstracts of the Forty-ninth Interscience Conference on Antimicrobial Agents and Chemotherapy, San Francisco, 2009*. Abstract H-899/378. American Society for Microbiology, Washington, DC, USA.